


VACUMSOL

Suction units for minor surgery

VACUMSOL AS60 – VACUMSOL AS65 – VACUMSOL AS65H

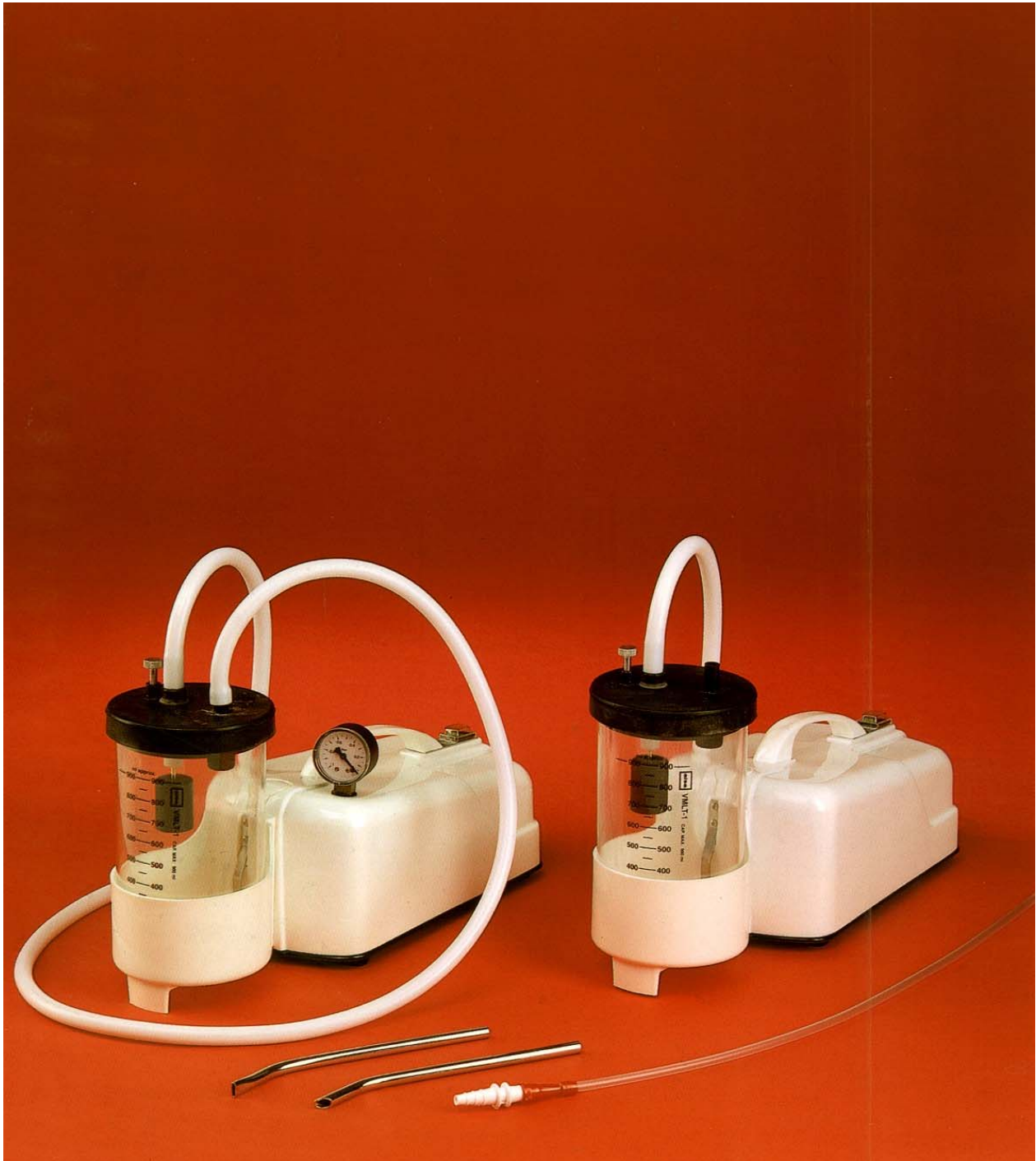

Quality System Certification



ISO 9001/EN 46001

alsa

CE 0051

Vacumsol models are three small portable suction units for minor surgery that differ one another in their working and suction power, and are therefore suitable for all possible uses: in home care after tracheotomy, hospital.

TECHNICAL FEATURES

- Cabinet made of shock proof plastic material
- Diaphragm aspirating group without maintenance (working noise less than 50 dB)
- Unbreakable and autoclavable bottle of lt. 1 made of polycarbonate with overflow device and vacuum regulator
- Manufactured according to the International Safety Standards IEC 601-1 and UNI EN ISO 10079/1
- Class II – type BF
- Weight: kg 2.7.

VACUMSOL AS/60

Maximum vacuum 650 mm/Hg (0.85 bar, 0.85 atm) displacement 17 l/min.

It is the base model of the series, with supply voltage, it is silent and easy to use and it is suitable for home care after tracheotomy, for neonatal suction...

VACUMSOL AS/65

Maximum vacuum 700 mm/Hg (0.9 bar, 0.9 atm) displacement 19 l/min.

It is more powerful than model VACUMSOL AS 60, it grants a better working and it is indicated also for pharyngeal suction, for ORL, for small surgery...

VACUMSOL AS/65-H

Technically identical to mod. AS 65, but equipped with vacuum-gauge.

It is the most suitable model for professional and hospital needs.

ACCESSORIES ON REQUEST

- Autoclavable bacterial filter (reusable and sterilizable)
- Saliva sucking cannula
- Mucus sucking cannula
- Soft disposable cannula CH14 (l. cm 60 Ø mm 4.7)
- Soft disposable cannula CH16 (l. cm 60 Ø mm 5.4)
- Soft disposable cannula CH18 (l. cm 60 Ø mm 6.1)
- Connection for soft cannulae.

ALSA

Apparecchi medicali s.r.l.

Via C. Bonazzi, 16

40013 CASTEL MAGGIORE (BO) ITALY

Tel. (+39) 051 700101 (r.a.) Fax 051 702182

Sito web: www.alsamed.com

e-mail: alsa@alsamed.com